
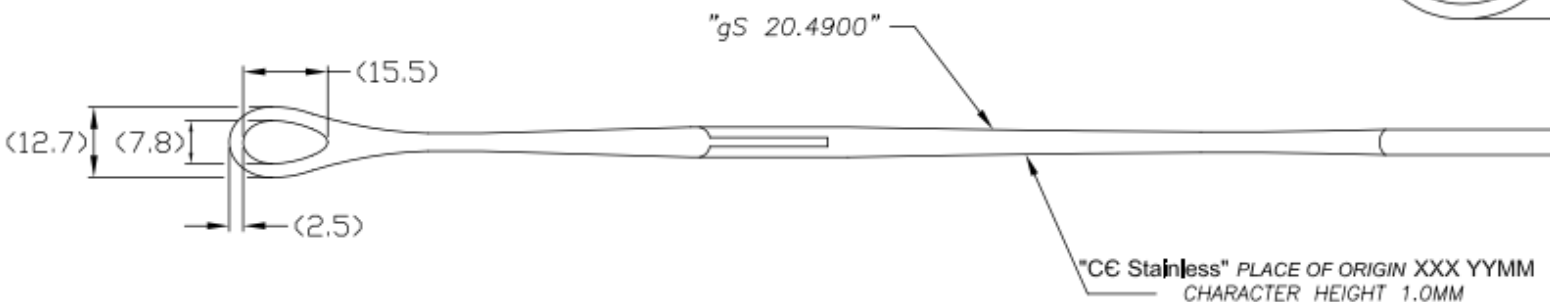
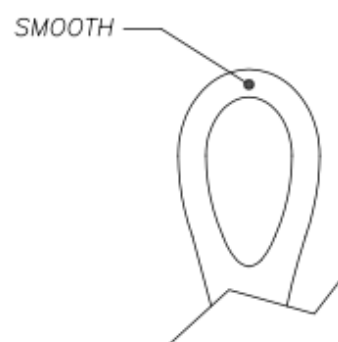


PART #	"A"
gS 20.4700	180
gS 20.4900	247

**FREIGEGERBEN
RELEASED**

Released 5.10.18

REVISIONS			
REV	DATE	APPROVED	DESCRIPTION
A	3/14/04	<i>JB</i>	NEW DRAWING
B	8/15/12	<i>JB</i>	UPDATE MARKINGS, NOTES AND PAPER SIZE; ADD FAMILY; SHOW SMOOTH DETAIL

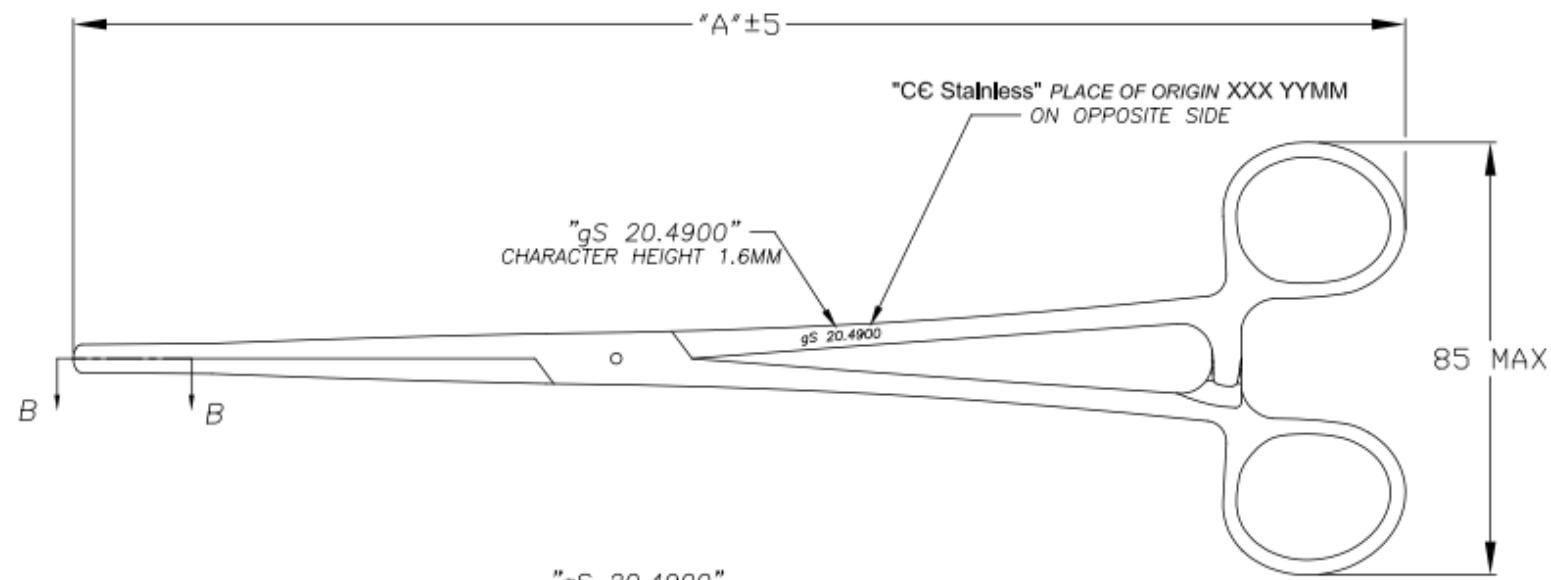

1. REMOVE ALL BURRS. NO SHARP EDGES UNLESS SPECIFIED.
2. PASSIVATE STAINLESS STEEL MATERIAL PER ASTM A967 OR EQUIVALENT.
3. MARKING TO BE ELECTRO CHEMICALLY ETCHED, FOLLOWED BY SPOT PASSIVATION.
4. MATERIAL HARDNESS: HRC 46±4.
5. DIMENSIONS SHOWN WITH INSTRUMENT ON FIRST TOOTH.
6. PART SHOWN IS 20.4900.

NOTES

THE INFORMATION CONTAINED HEREIN IS CONSIDERED PROPERTY OF gSource AND MAY NOT BE DISCLOSED FOR ANY PURPOSE OTHER THAN gSource RELATED MANUFACTURING REQUIREMENTS WITHOUT WRITTEN CONSENT FROM gSource.

gSource			
DO NOT SCALE			
TOLERANCE UNLESS SPECIFIED OTHERWISE X. ± 1 X.X ± 0.5 X.XX ± 0.12 ANGLES ± 2°	METRIC SCALE: 1:1	THIRD ANGLE PROJECTION SIZE: A3 SHEET: 1 of 1	
DRAWN: <i>Jim B.</i>	APPROVED: <i>JB</i>	DATE: 3/14/04	
MATERIAL: 1.4021		FINISH: SATIN	
DESCRIPTION: FOERSTER FORCEPS STRAIGHT, SMOOTH			
DWG/PART# 20.4900			REVISION: B